## INSTRUCTIONS

EN

Intended for use in the general population.

English

1. **Select** the area for the placement of the Grip-Lok. Note: catheter, tube or line should be already placed/inserted.
2. **Prepare** the skin according to the standard hospital protocol for dressing application.
3. **Open** the top flap and slide the Grip-Lok under the catheter, tube or line, centering it above the adhesive area.
4. **Pull** the paper backing from one side of the Grip-Lok, then the other, to secure in desired position on the skin.
5. **Remove** the interior liner to expose the adhesive areas.
6. **Place** the catheter, tube or line in the center of the exposed adhesive.
7. **Secure** top flap over the catheter, tube or line and apply gentle pressure to the center fabric area over the adhesive.

**Direction:** for added skin protection, use skin prep pad prior to application of stabilization device.

**Note:** Use of an alcohol swab may help removal of securement device from the skin.

⚠ If you experience severe redness, itching, swelling, or irritation of the skin, consult your physician as this may be a sign of an allergic reaction.

⚠ Re-use of this device may change its mechanical or biological features and may cause device failure, allergic reactions or bacterial infections.

⚠ Replace securement device if soiled or saturated in fluid or if device shows signs of wear or damage.

**Note:** If being used in a medical facility, dispose according to medical facility and municipal guidelines.

⊘ Do not use if the product's sterile barrier system or its packaging is compromised.

Report any serious incident related to device use to TIDI Products and the Member State competent authority.

HYPALLERGENIC

TIDI Products, LLC  
570 Enterprise Drive  
Neeenah, WI 54956, USA  
Phone: + 1.920.751.4300  
tidiproducts.com/securement

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UK  
CA  
0086

UK Responsible Person:  
Emergo Consulting (UK) Limited  
c/o Cr 360 - UL International  
Compass House, Vision Park Histon  
Cambridge CB24 9BZ  
United Kingdom

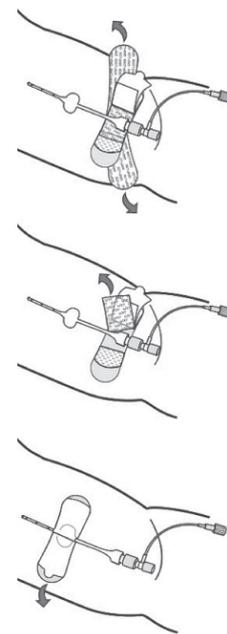
Grip-Lok is a registered trademark of TIDI Products, LLC.

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Grip-Lok är ett registrerat varumärke som tillhör TIDI Products, LLC.
Grip-Lok es una marca registrada de TIDI Products, LLC.