3.11"
This panel visible when folded

3.11"
This panel visible when folded

3.11"

Posey® 8269 Alarm Attachment Bracket Application Instructions

NAME:

REF 8269 Alarm Attachment Bracket

DESCRIPTION OF PRODUCT:
Posey Alarm Mounting Bracket with Attachment Strap and Quick Release Buckle


INTENDED USE: The 8269 Alarm Attachment Bracket is suitable for bed and chair/wheelchair mounting.

CONTRAINDICATIONS: None

INSTRUCTIONS:

1. Choose location on headboard/footboard or chair/wheelchair where there is an opening and the strap can be attached so that the patient cannot reach or tamper with the alarm or connections.
2. Buckle strap onto bed (Fig. 1) or chair/wheelchair and tighten strap securely.
3. Slide alarm onto bracket from top down until it is firmly in place (Fig. 2). Make sure alarm indictor lights are in clear view of staff.
4. To remove alarm, gently push release lever IN while sliding alarm up and out (Fig. 3).

DISPOSAL: Dispose of the bracket in accordance with the applicable local or facility guidelines. Report any serious incident related to device use to TIDI Products and the Member State competent authority.

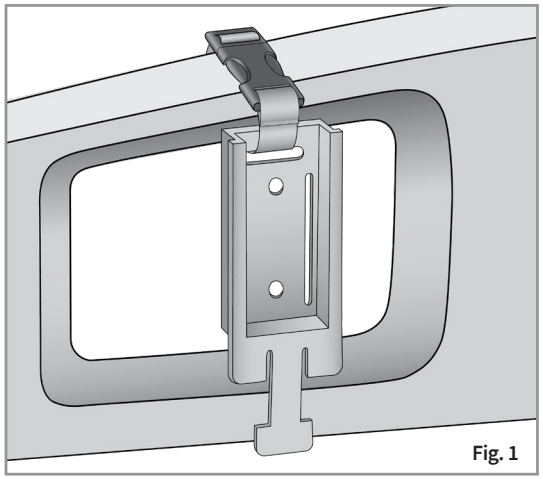
Fig. 1

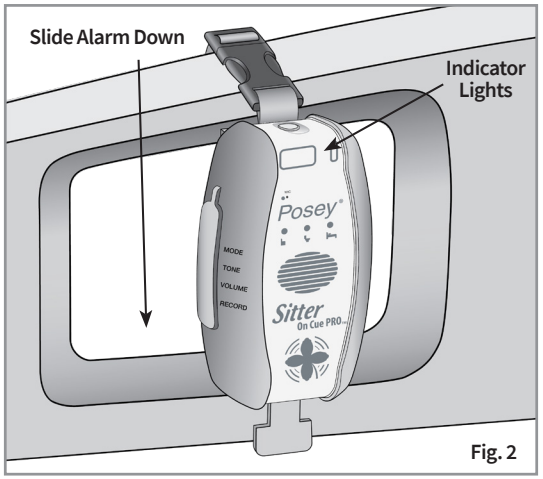


Fig. 2

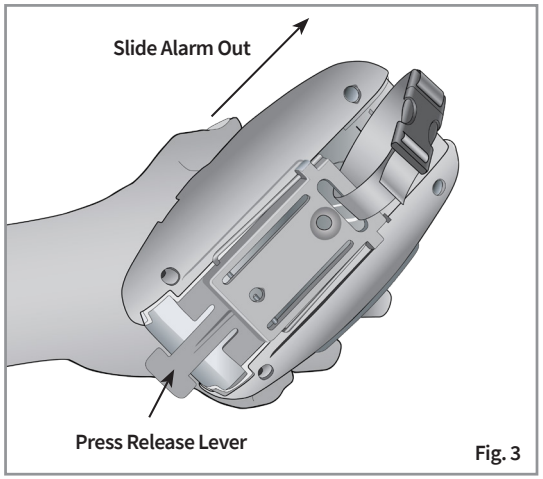


Fig. 3

TIDI Products, LLC • 570 Enterprise Drive, Neenah, WI 54956 USA
Phone: 1.800.521.1314 • International: +1.920.751.4036 • www.tidiproducts.com

© 2026 TIDI Products, LLC. All rights reserved. Posey is a registered trademark of TIDI Products, LLC.

TIDI Support Caregivers. Protect Patients.

M7001 REV C 2026-04-02

Printer: please note dashed fold lines. The upper left and upper middle sections should be visible when document is fully folded.
Folded size will be ~3.11 x 3.375