
Posey® 8208CS Command™ Strip Alarm Mounting System Application Instructions

NAME:

REF 8208CS Command™ Strip Alarm Mounting System

DESCRIPTION OF PRODUCT:
Posey Command™ Strip Alarm Wall Mounting Bracket

INTENDED USE:
To mount Posey Alarms to hospital room wall without the need to drill holes

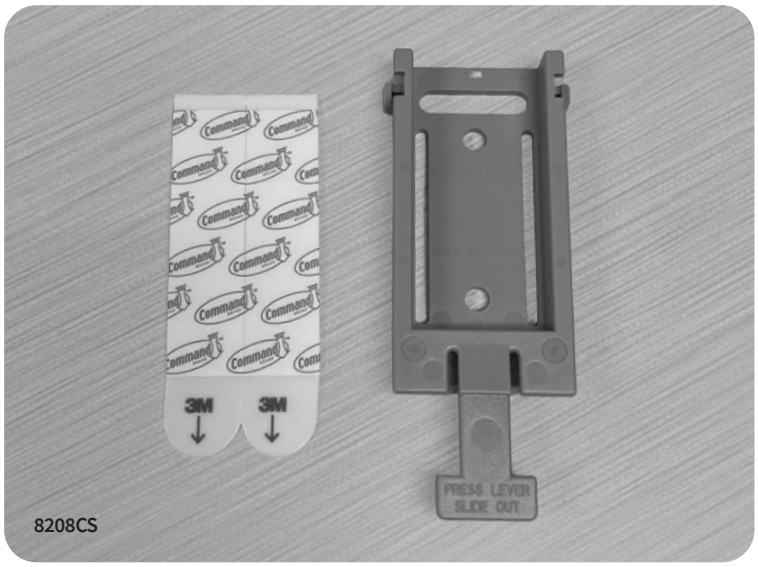

CONTRAINDICATIONS: None


INSTRUCTIONS: Installing the 8208CS to Wall Surface

1. Use isopropyl alcohol and a clean cloth to clean the back of the Posey wall bracket and the wall surface where bracket will be attached. Wipe gently and let dry.

⚠ WARNING DO NOT USE any household or industrial cleansers as they may leave a residue which will adversely affect adhesion.




Isopropyl Alcohol

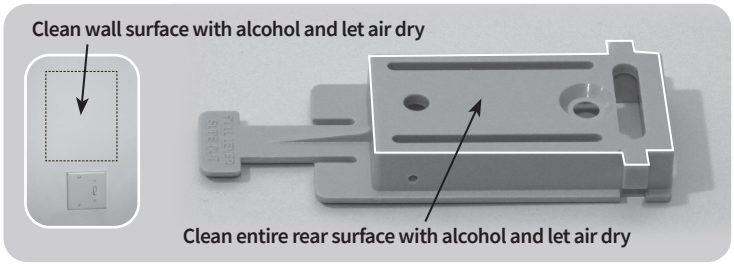



2. Each Command™ Strip has a "RED" release liner and "BLACK"release liner on opposite sides of the strip (Fig. A).

3. Remove RED release liner from Command™ Strips.

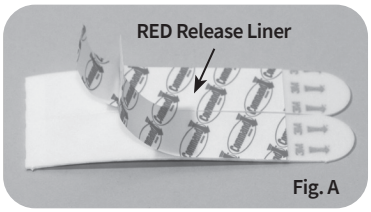


4. Align the adhesive with the top edge of wall bracket and press adhesive firmly in place against the bracket. Apply pressure for 10 seconds to activate the pressure sensitive adhesive (Fig. B).

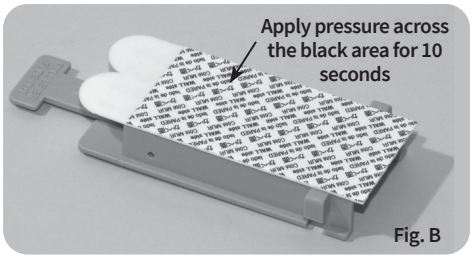

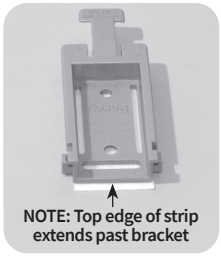

5. Remove BLACK release liners and secure bracket to wall. Press bracket against wall for 30 seconds to activate the pressure sensitive adhesive. Wait one (1) hour before use (Fig. C).



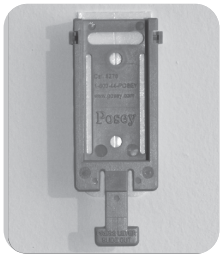

⚠ WARNING The Posey 8208CS will adhere to the following types of wall surfaces: painted, stained or varnished wood, glass tile, painted cinder block, plaster, metal or wallboard. It is not suitable for wallpaper or brick. Surface temperatures should be above 50°F (10°C) for proper adhesion and adhesion could lose strength at temperatures above 105° (40°C).


Inspect all components for damage prior to use.

REMOVAL INSTRUCTIONS:

1. Locate the two Command™ Strip release tabs near the bottom of the bracket. If needed, use a non-sharp object to help space them off the wall to aid in grabbing them.

2. Hold the wall bracket with one hand and pull the release tabs with your other hand.



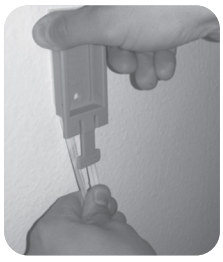



⚠ WARNING DO NOT pull release tabs towards you. Pull straight down along the wall at least 15 inches. The strips will stretch and eventually release from the wall.

DISPOSAL: Dispose of the bracket in accordance with the applicable local or facility guidelines. Report any serious incident related to device use to TIDI Products and the Member State competent authority.

TIDI Products, LLC • 570 Enterprise Drive, Neenah, WI 54956 USA
Phone: 1.800.521.1314 • International: +1.920.751.4036 • www.tidiproducts.com

© 2026 TIDI Products, LLC. All rights reserved. Posey is a registered trademark of TIDI Products, LLC. 3M and Command™ are the worldwide trademarks or registered trademarks of 3M.

M6162 REV C 2026-03-19