
Posey® 8208L Alarm Locking Bracket Application Instructions

NAME:

REF 8208L Posey Alarm Locking Bracket

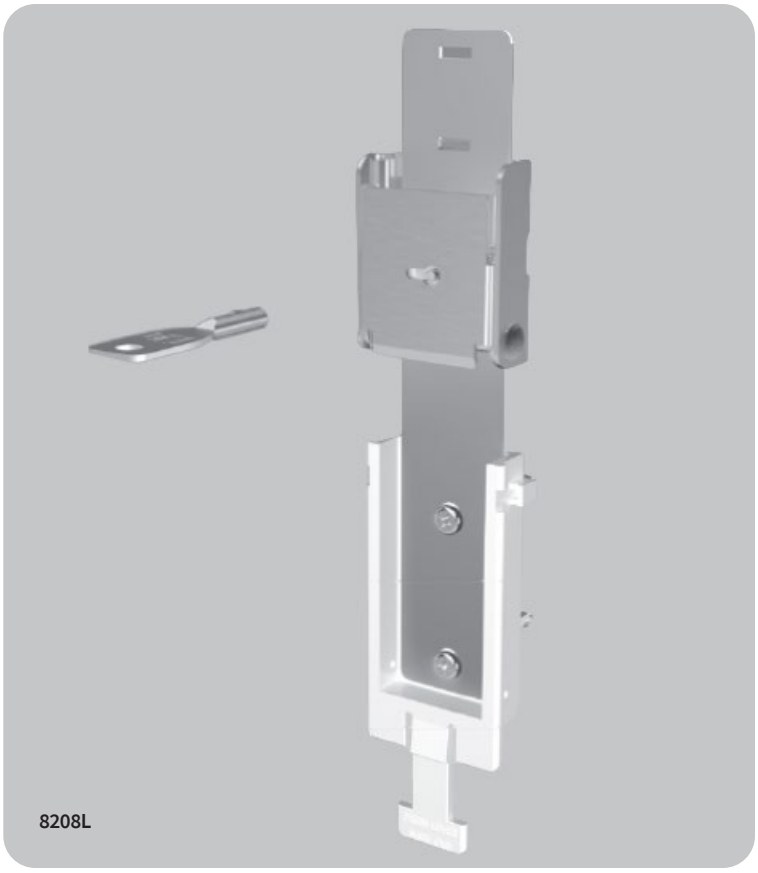

DESCRIPTION OF PRODUCT: Posey Alarm Locking Bracket

INTENDED USE: To mount and lock Posey alarms to wall bracket 8208

CONTRAINDICATIONS: None

INSTRUCTIONS:


⚠ WARNING The 8208L is designed to mount to a flat wall surface only.


1. Choose a wall location out of the patient's reach, where the alarm's indicator lights will be in clear view of staff.
2. Screw attachment:
   - Position wall bracket with back (flat side) against wall
   - Using bracket as a template, mark spots to insert the screw anchors into the wall. (Fig. 1)

   ⚠ WARNING Make sure it is safe to drill and there are no pipes or electrical wires that could be damaged.

   ⚠ WARNING If using your own mounting hardware, use a screw with a head diameter (or washer) large enough to secure wall bracket/Narrow Vertical Slide to wall.
3. Drill holes where marked and install anchors.
4. Place the Narrow Vertical Slide into the Wall Bracket so that the bottom of the Slide is flush against the bottom wall of the bracket. (Fig. 2a)
5. Position Narrow Vertical Slide and Wall Bracket assembly over wall anchor holes; insert and tighten screws. (Fig. 2b)
6. Slide alarm onto bracket from top down until it is firmly in place. (Fig. 3a)
7. Slide buckle lock on Narrow Vertical Slide, to top of alarm. (Fig. 3b)
8. Lock buckle by pressing buckle lock to engage the slot in Narrow Vertical Slide. (Fig. 4)
9. To remove alarm, first unlock and remove buckle lock. (Fig 5a-b) Then gently push bracket release lever IN (Fig 5c) while sliding alarm up and out. (Fig 5d)
10. Replace buckle lock and lock onto Narrow Vertical Slide when not in use.

DISPOSAL: Dispose of the bracket in accordance with the applicable local or facility guidelines. Report any serious incident related to device use to TIDI Products and the Member State competent authority.

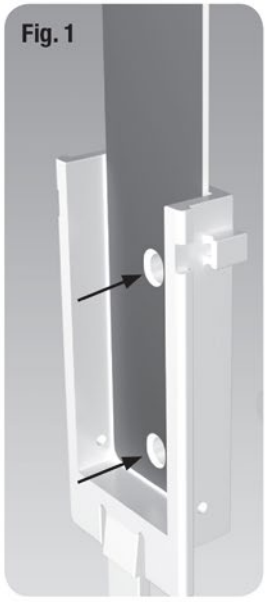

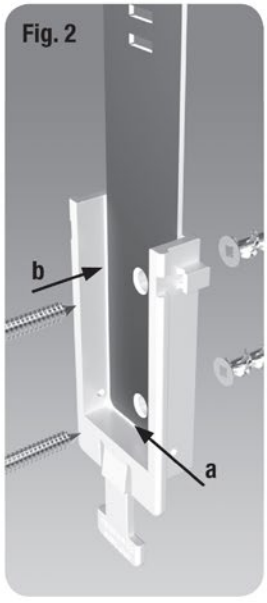







TIDI Products, LLC • 570 Enterprise Drive, Neenah, WI 54956 USA
Phone: 1.800.521.1314 • International: +1.920.751.4036 • www.tidiproducts.com

© 2026 TIDI Products, LLC. All rights reserved. Posey is a registered trademark of TIDI Products, LLC.

TIDI Support Caregivers. Protect Patients.

M6163 REV D 2026-03-19